





It's a fact. Clinton Pediatric Equipment is among the best in the medical equipment business for innovation, design superiority, breadth of line and durable performance. And surprisingly enough, it is among the most affordable.



Our *Clinton Kids™* products feature exclusive graphics. Designed by Clinton illustrators, these colorful graphics create a comfortable, kid-friendly environment for children of all ages.



Completely coordinated. Completely cost effective. Completely easy. Clinton has done the design work for you with our total designed exam room solutions — *Clinton Complete*. Clinton Complete provides package pricing on the purchase of specific tables and coordinating cabinets. Plus, we recommend paint colors, flooring and floor plan configurations. Add coordinating graphics and accessories to achieve a total *Clinton Complete* look.

Look for *Clinton Complete* designed solutions throughout this catalog.

We invite you to examine our products. Imagine the possibilities, and select the equipment, from the one source that is committed to your pediatric needs.







## FUN SERIES SCALE TABLES

The Clinton Kids™ Fun Series Scale Tables are the perfect addition to any pediatric office or clinic. Fun shapes are combined with an accurate digital scale and soft safety top for the ultimate in performance and patient comfort.

### Fun Series Scale Tables feature:

- Comfortable work height of 36" (91.44 cm)
- Top rimmed on 3 sides for safety
- Beveled edges & rounded corners on top
- Easy-clean, replaceable center pad
- Built-in infantometer
- Easy access to weight measurement and storage
- Built-in concealed paper dispenser & cutter
- Easy-clean laminate surfaces
- Ample 2-door storage
- Optional storage locks available (-055)

### Doran® Scales feature:

- Dependable, accurate Doran® scale
- Easy-read digital display
- Auto calibration on start-up
- Auto shut-off
- Electric or battery operation
- English or Metric weights
- Maximum weight 40 lbs. (18.14 kg)



top rimmed on 3 sides



beveled sides & rounded corners



infantometer



digital readout

### Zoo Bus Scale Table

- 1 end door and 1 side door storage with 1 adjustable shelf
- Graphics on 1 side, front and back

7822	top length	top width	
	43"	24"	
	109.22 cm	60.96 cm	
	base length	base width	base height
	47"	21"	36"
	119.38 cm	53.34 cm	91.44 cm



back detail



ZOO BUS SCALE COMPLETE includes 7822 Zoo Bus Scale Table & 6120-BW Coordinating Cabinets  
ADD:  
COORDINATING GRAPHICS shown here for a total "Zoo Bus" exam room designed environment

For special package discount pricing, ORDER USING MODEL # 7822-X



**Engine K-9 Scale Table**

- 1 end door and 1 side door storage with 1 adjustable shelf
- Graphics on 1 side, front and back

<b>7833</b>	<i>top length</i>	<i>top width</i>	
	43"	24"	
	109.22 cm	60.96 cm	
	<i>base length</i>	<i>base width</i>	<i>base height</i>
	47"	21"	36"
	119.38 cm	53.34 cm	91.44 cm



back detail

**Outback Buggy Scale Table**

- 1 end door and 1 side door storage with 1 adjustable shelf
- Graphics on 1 side, front and back

<b>7844</b>	<i>top length</i>	<i>top width</i>	
	43"	24"	
	109.22 cm	60.96 cm	
	<i>base length</i>	<i>base width</i>	<i>base height</i>
	47"	21"	36"
	119.38 cm	53.34 cm	91.44 cm



**OUTBACK BUGGY SCALE COMPLETE**  
–discount package pricing when  
combining 7844 table with coordinating  
cabinets

**7844-X** Must order using this model # for package discount



back detail



ENGINE K-9 SCALE COMPLETE includes  
7833 Engine K-9 Scale Table &  
6130-BW Coordinating Cabinets  
ADD:  
COORDINATING GRAPHICS shown  
here for a total "Engine K-9"  
exam room designed environment

For special package discount pricing, ORDER USING MODEL # 7833-X



## SELECT SERIES SCALE TABLES

The Clinton Select Series Scale Tables offer a wide range of features, and are available in over a dozen laminate colors for the front. Add pediatric fun by customizing with exclusive Clinton Select Series Graphics (see page 63).

Select Series Scale Tables feature:

- Comfortable work height of 36" (91.44 cm)
- Top is rimmed on three sides for safety
- Top's front and back rim edges are beveled and the corners rounded for safety
- Easy-clean, replaceable center pad
- Built-in infantometer with folding ends for measuring height
- Easy access to weight, measurement and storage
- Built-in concealed paper dispenser & cutter
- Recessed front toe kick
- Adjustable built-in levelers provide a stable surface
- Easy-clean all laminate base
- Ample storage with various door and drawer configurations
- Optional storage area locks available

See page 4 for Doran® scale features.

### Pediatric Scale/Treatment Table

- Features 2 drawers and 2 doors

7840	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm



### Pediatric Scale/Treatment Table

- Features 1 drawer and 2 open shelves

7810	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

### Pediatric Scale/Treatment Table

- Features 2 doors with 1 adjustable shelf


7820	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm



### Pediatric Scale/Treatment Table

- Features 2 doors and 2 drawers

7830	length	width
	43"	24"
	109.22 cm	60.96 cm

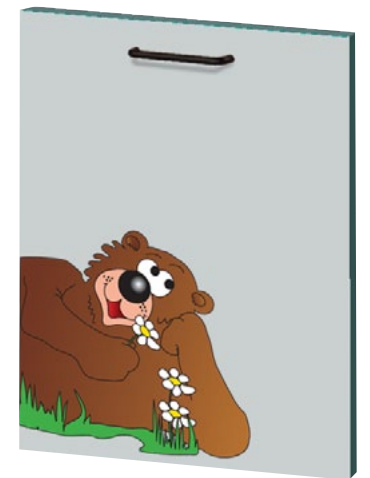

Add Select Series Graphics to doors and drawer fronts for a fun custom look (see page 63)





## IMAGINATION SERIES SCALE TABLES

The Clinton Kids™ Imagination Series Scale Tables take office decor to a new dimension with bright illustrations and the ultimate in comfort, performance and accuracy.

*Imagination Series Scale Tables feature:*

- Comfortable work height of 36" (91.44 cm)
- Top rimmed on 3 sides for safety
- Beveled edges & rounded corners on top
- Easy-clean, replaceable center pad
- Built-in infantometer
- Easy access to weight measurement and storage
- Built-in concealed paper dispenser & cutter
- Easy-clean laminate surfaces
- Adjustable built-in levelers provide stability
- Ample 2-door and 1 drawer storage
- Optional storage locks available (-055)
- All models feature **Clinton Complete™** pricing when combined with coordinating cabinets

See page 4 for Doran® scale features.

### Dino Days Scale Table


- 2-door and 1 drawer storage

7834	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm



right side detail

left side detail

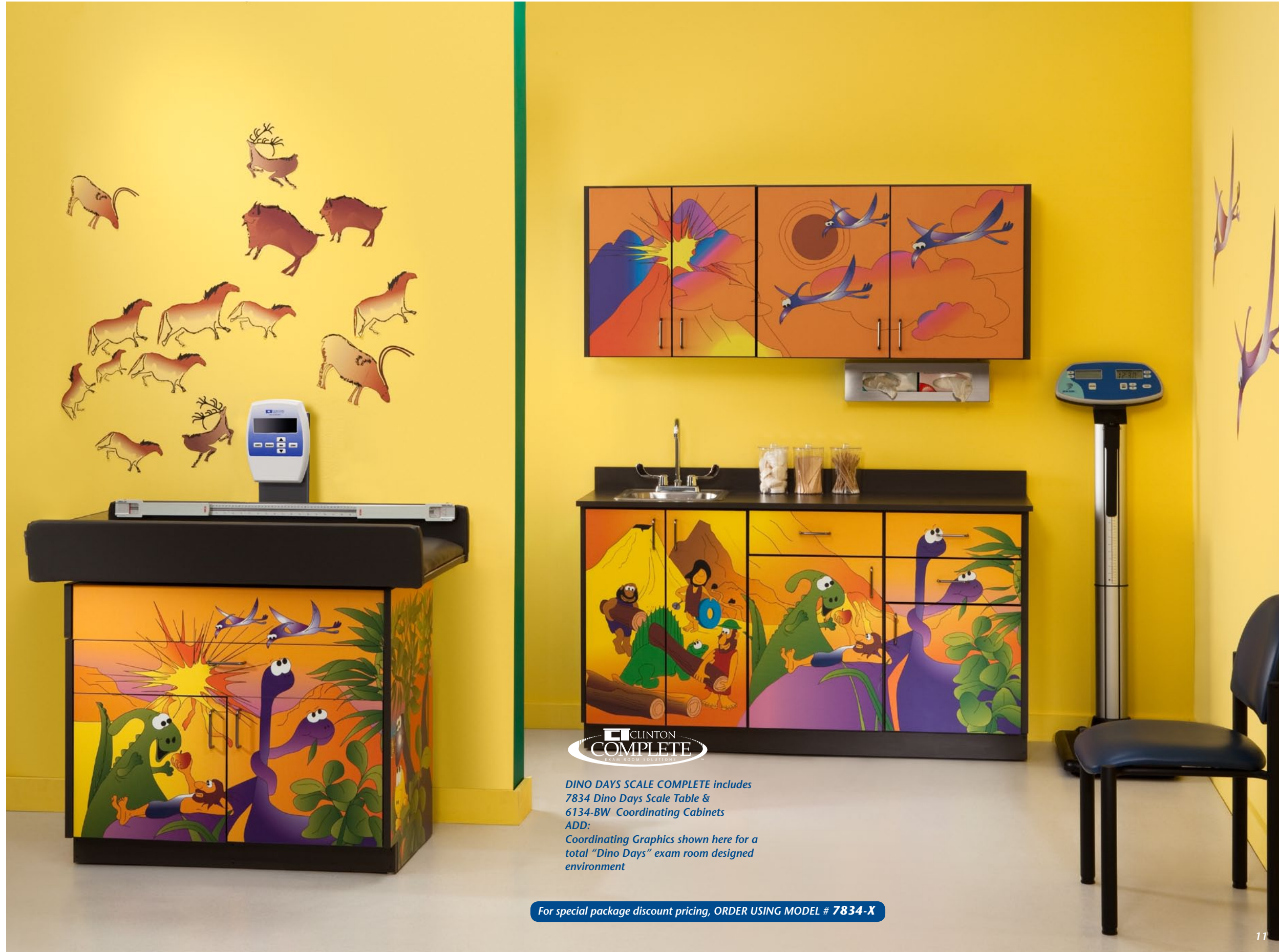

**DINO DAYS SCALE COMPLETE** includes  
7834 Dino Days Scale Table &  
6134-BW Coordinating Cabinets  
ADD:  
Coordinating Graphics shown here for a  
total "Dino Days" exam room designed  
environment

For special package discount pricing, ORDER USING MODEL # 7834-X



**Cool Pals Scale Table**

- 2-door and 1 drawer storage

7831	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

*COOL PALS SCALE COMPLETE—discount package pricing when combining 7831 table with coordinating cabinets (see page 44 for cabinets and features)*

**7831-X** Must order using this model # for package discount



right side detail

left side detail



**Rainforest Follies Scale Table**

- 2-door and 1 drawer storage

7832	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

*RAINFREST FOLLIES SCALE COMPLETE—discount package pricing when combining 7832 table with coordinating cabinets (see page 44 for cabinets and features)*

**7832-X** Must order using this model # for package discount



right side detail

left side detail



**Space Place Scale Table**

- 2-door and 1 drawer storage

7835	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

*SPACE PLACE COMPLETE—discount package pricing when combining 7835 table with coordinating cabinets (see page 45 for cabinets and features)*

**7835-X** Must order using this model # for package discount



right side detail

left side detail



**Ocean Commotion Scale Table**

- 2-door and 1 drawer storage

7836	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

*OCEAN COMMOTION SCALE COMPLETE—discount package pricing when combining 7836 table with coordinating cabinets (see page 42 for cabinets and features)*

**7836-X** Must order using this model # for package discount



right side detail

left side detail





**Alley Cats & Dogs Scale Table**

- 2-door and 1 drawer storage

7837	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm



right side detail



left side detail



**ALLEY CATS & DOGS SCALE COMPLETE**  
includes 7837 Alley Cats & Dogs Scale  
Table & 6137-BW Coordinating Cabinets  
ADD:  
Coordinating Graphics shown here for  
a total "Alley Cats & Dogs" exam room  
designed environment



For special package discount pricing, ORDER USING MODEL # 7837-X



**Serengeti Sunrise Scale Table**

- 2-door and 1 drawer storage

7851	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

**SERENGETI SUNRISE SCALE COMPLETE**—discount package pricing when combining 7851 table with coordinating cabinets (see page 46 for cabinets and features)

**7851-X** Must order using this model # for package discount



right side detail

left side detail



**Kangaroo Country Scale Table**

- 2-door and 1 drawer storage

7853	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

**KANGAROO COUNTRY SCALE COMPLETE**—discount package pricing when combining 7853 table with coordinating cabinets (see page 47 for cabinets and features)

**7853-X** Must order using this model # for package discount



right side detail

left side detail



**Southwestfest Scale Table**

- 2-door and 1 drawer storage

7852	length	width	height
	43"	24"	36"
	109.22 cm	60.96 cm	91.44 cm

**SOUTHWESTFEST SCALE COMPLETE**—discount package pricing when combining 7852 table with coordinating cabinets (see page 47 for cabinets and features)

**7852-X** Must order using this model # for package discount



right side detail



left side detail





## WEIGHING PROCESS



### Weighing a Patient

1. Place the patient on the scale platform.
2. The indicator will automatically turn on and display the weight.
3. The indicator will beep and hold the weight on the display until the scale automatically shuts off.

**NOTE:** If the scale does not automatically turn on, remove the patient and press zero.

### IMPORTANT:

- The entire exam table top is part of the scale. Leaning against or touching the table top while weighing a patient will result in weight errors. If the top was compromised during weighing, press the REWEIGH button to get another readout.
- If the patient on the table is moving and cannot remain reasonably still, an inaccurate readout may occur, press the REWEIGH button to weigh the patient again.
- When using the storage compartment located on the top left side of the cabinet, take care that items are not in contact with underside of the top. This is part of the scale and may result in inaccurate readout.

## MEASURING PROCESS



### Using the Infantometer

Clinton Scale tables are equipped with a handy infantometer that is also very simple to use. The ruler has centimeters on top and inches on the bottom for your convenience.

1. If not already locked into position, lock the left paddle to the left edge of the ruler. The left end factors in the length of the slider so the protruding paddle is always zero.
2. Gently place the infant between the paddles.
3. Gently slide the right paddle along the ruler to the infant's feet.
4. The red arrow on the right paddle indicates the accurate length of the infant, factoring in the width of the slider.





*Ask about Clinton's entire line of products.  
Call for any of our product catalogs today.*

- Medical Tables
- Medical Casework
  - Power Tables
- Seating & Medical Accessories
  - Phlebotomy Equipment
  - Physical Therapy Equipment



800-441-9131 • 717-848-3519

Toll-Free Fax 866-522-1233

Fax 717-843-5871

[www.clinton-ind.com](http://www.clinton-ind.com)

Scan with your smartphone  
to see our Pediatric Equipment video

