

# AV4028

Gynaecologic bed, 3 sections, 2 articulations, fixed height

Gynaecological couch, fixed height, very handy thanks to its ergonomics.  
Adjustable backrest with semiautomatic racks, maximum slope 70°, fixed pelvis section.  
Removable foot rest with 2 positions: vertical or horizontal; the central hollow allows the extraction of the basin placed under the seat section.  
The structure is composed by legs made of chromed tubular steel, connected to the main frame through melting of alloy.  
The main frame support the lying area, the couple of racks, the basin support and the clamps for the thighs supports.  
The lying area are made of laminated panels, padded and upholstered of plastic, waterproof fabric.  
The seals of the upholstery are double against tears.

### Equipped with


- stainless steel basin, extractable and detachable
- adjustable ABS thighs supports, fixed through lateral clamps, adjustable in height



### Accessori disponibili / Available accessories



**AV4088**  
Couple of lateral profiles for coupling of all accessories

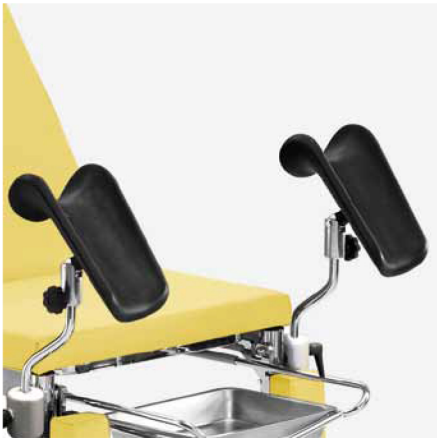

**AV4081**  
Thighs support made of soft polyurethane with steel inner support



**AV4082**  
Covering made of fireproof plastic class 1 IM (additional covering for AV4081)



**AV4098**  
I. V. stand with two hooks made of chromed steel with clamp for height adjustment, to be hooked on lateral profiles



**AV4097B**  
Paper roll support mounted on the bed base, with support made of chromed steel



**AV4074**  
Adjustable handles made of soft integral polyurethane, usable as small siderails

## Dati tecnici / Technical features

DATI DIMENSIONALI / DIMENSIONS	
Dimensioni totali / Total dimensions	185x60x80h cm
Dimensioni piano di coricamento / Lying area dimensions	185x60 cm
Dimensioni schienale / Back rest dimensions	90x60 cm
Dimensione bacino / Basin section dimensions	42x60 cm
Dimensione gambale / Foot rest dimensions	54x60 cm
Inclinazione schienale / Back rest tilt	70°
Portata / Capacity	150 kg
Peso / Weight	34 kg
Volume / Volume	0,260 mc
Spessore supporto imbottitura / Padding support thickness	18 mm
Densità imbottitura / Padding density	30 kg/mc
Rivestimento ignifugo classe / Fireproof covering class	1 IM
Dispositivo medico marcato CE DIR.93/42 appartenente alla classe 1 / Medical device marked CE DIR. 93/42 belonging to class 1	