

AV4030

Couch for medical centre

Three sections couch, to be adapted for gynaecology as well.
Transformed simply by lifting the backrest or the central section through a six positions rack.
Foot/pelvis section is fixed, the central hollow allows the extraction of the basin placed under the seat section; a kind of “closing section” closes the hollow when the basin is not used.
The structure is composed by legs made of chromed tubular steel, connected to the main frame through melting of alloy.
The main frame supports the lying area, the couple of racks, the basin support and the clamps for the thighs supports as well.
The three lying areas are made of laminated panels, padded and upholstered of plastic, waterproof fabric.
The seals of the upholstery are double against tears.

Equipped with


- stainless steel basin, extractable and detachable
- adjustable ABS thighs supports, fixed through lateral clamps, adjustable in height



Accessori disponibili / Available accessories



AV4088
Couple of lateral profiles for coupling of all accessories

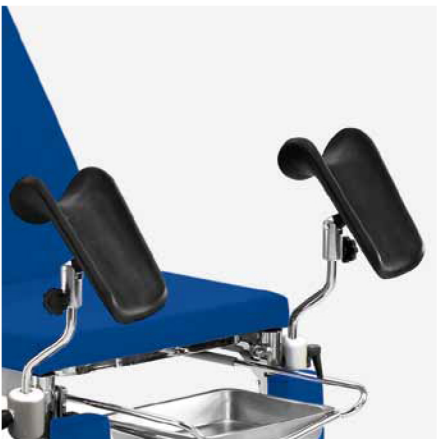

AV4081
Thighs support made of soft polyurethane with steel inner support



AV4082
Covering made of fireproof plastic class 1 IM (additional covering for AV4081)



AV4098
I. V. stand with two hooks made of chromed steel with clamp for height adjustment, to be hooked on lateral profiles



AV4097B
Paper roll support mounted on the bed base, with support made of chromed steel



AV4074
Adjustable handles made of soft integral polyurethane, usable as small siderails.

Dati tecnici / Technical features

DATI DIMENSIONALI / DIMENSIONS	
Dimensioni totali / Total dimensions	186x70x80h cm
Dimensioni piano di coricamento / Lying area dimensions	186x70 cm
Dimensioni schienale / Back rest dimensions	62x70 cm
Dimensioni sezione centrale / Central section	85x70 cm
Dimensione pediera-bacino / Footrest-pelvis section	39x70 cm
Inclinazione schienale / Back rest tilt	70°
Portata / Capacity	150 kg
Peso / Weight	35 kg
Volume / Volume	0,3 mc
Spessore supporto imbottitura / Padding support thickness	18 mm
Densità imbottitura / Padding density	30 kg/mc
Rivestimento ignifugo classe / Fireproof covering class	1 IM
Dispositivo medico marcato CE DIR.93/42 appartenente alla classe 1 / Medical device marked CE DIR. 93/42 belonging to class 1	