

Serie BS altezza fissa

BS Line - fixed height


Base composed by two removable sections, made of steelectro galvanized, hooked on plastic supports.
At the four corners bumper bars in soft plastic material, non marking, andholes for push handles and accessories.

At the head side supports for drainage bags andpredisposition for oxygen bottle holder.Backrest is moved by bilateral lever and assisted in the movement withpneumatic spring.

Lightweight frame made of thick steel profiles andepoxy paint.
Four castors Ø 200 mm, one swivel, one fixed directional ;two swivel castors at head side with spring brake always on (BS1602) these can be unlocked operating on the handle on the head side.

The device has been designed and built by GIVAS SRL in accordance to class I safety requirements provided by 93/42/CEE Legislation of 14th June 1993 Council, and further mod. integrations, about medical devices.

Versione dotata di freno costante
Model with constant brake

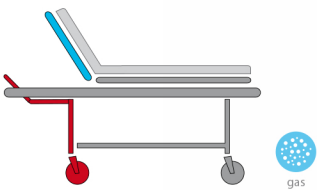

MOVIMENTI / MOTIONS



Kit: BSK0005

La barella raffigurata è composta dagli articoli:
The depicted stretcher is composed of the following items:

- **BS1600** - base barella a 2 sezioni ad altezza fissa
2-section stretcher (frame only) with fixed height
- **BZ1604** - coppia sponde a ribalta
folding side rails
- **BZ1501** - materassino rivestito antistatico
antistatic mattress
- **ZBZ15241** - asta flebo a 4 ganci regolabile in altezza
adjustable IV stand with 4 hooks

- **BZ1520** - portabombola in acciaio cromato
dim. int. Ø 155x410 mm
oxygen bottle holder made of chromed steel,
inner measures Ø 155 x 410 mm
- **BZ1525** - supporto reggirotolo in acciaio inox
per lenzuolino
Paper support made of stainless steel
for sheet
- **BZ1606** - cestello porta indumenti
clothes basket
- **BZ1535** - N° 2 maniglioni di spinta con sedi asta flebo
e portamonitor
push handles with holes for I.V. stand
and shelf for monitor



EQUIPAGGIAMENTO | EQUIPMENT

COMPOSIZIONE | COMPOSITION

BS1602	2 sezioni, altezza fissa, freno costante 2 sections, fixed height, constant brake	■
	Ruote antistatiche Ø 200 mm Antistatic castors Ø 200 mm	■
	Paracolpi Bumpers	■
BS1600	2 sezioni, altezza fissa, freno indipendente 2 sections, fixed height, independent brake	■
	Ruote antistatiche Ø 200 mm Antistatic castors Ø 200 mm	■
	Paracolpi Bumpers	■

COMPONENTI OPZIONALI | OPTIONAL COMPONENTS

BZ1604	Coppia sponde a ribalta Folding side rails	●
BZ1607	Coppia sponde a compasso Folding side rails with compass movement	●
BZ1501	Materassino antistatico Antistatic mattress (per/for BZ1604)	○
BZ1506	Materassino antistatico Antistatic mattress (per/for BZ1607)	○
BZ1535	Maniglione di spinta Push handle	○
BZ1518	Portamonitor Monitor support	○
BZ1524	Asta flebo I.V. stand	○
BZ1520	Portabombola O ₂ Oxygen tank holder	○
BZ1525	Supporto per reggirotolo Paper holder	○
BZ1565	Piano radiotrasparente Radio transparent surface	●

DATI TECNICI DIMENSIONALI | DIMENSIONAL TECHNICAL DATA

Lunghezza totale Total length	mm	2040
Larghezza totale Total width	mm	700
Superficie di appoggio Lying surface	mm	1980x580
Altezza Height	mm	780
Inclinazione schienale Backrest incline	deg	80°
Spazio tra basamento e piano Space between basement and lying surface	mm	430
Materasso Mattress	cm	193x58x5
Peso totale Total weight	kg	48
Portata in sicurezza Safety load	kg	170

LEGENDA | KEY

- di serie | standard
- obbligatorio scelta | compulsory by choice
- opzionale assemblato | optional assembled
- opzionale collocabile | optional placed